

Cytokeratin 5/14

Mouse Monoclonal Antibody

Mob433

PDM140

PDM140-HL

****This data sheet is applicable to all sizes (volume) of product. Actual volume is indicated on vial.**

Document #: DS-0169-E
Effective Date: 12/18/2023

Immunogen	Clone	Species	Isotype	Primary Antibody Diluent
XM26: prokaryotic recombinant protein corresponding to 103 amino acid portion of C-terminal region of the human cytokeratin 5 molecule. LL002: Synthetic peptide of the extreme C-terminal of human keratin-14.	XM26 + LL002	Mouse	IgG1 & IgG3	K004

Lot specific Ig concentration available upon request.

Catalog #	Description
Mob433	concentrated antibody for use with Diagnostic BioSystems PolyVue™ Plus - Two Step Detection System
PDM140	ready to use antibody for use with Diagnostic BioSystems PolyVue™ Plus - Two Step Detection System
PDM140-HL	Ready to use antibody in RFID tagged vials for use on the HighLighter Staining System

Intended Use

For In Vitro Diagnostic Use. This product is intended for qualitative immunohistochemistry with normal and neoplastic formalin-fixed, paraffin-embedded tissue sections, to be viewed by light microscopy. Clinical interpretation of staining results should be accompanied by histological studies with proper controls. Patients' clinical histories and other relevant diagnostic tests should be utilized by a qualified person(s) when evaluating and interpreting results.

Summary and Explanation

This antibody cocktail reacts with cytokeratin 5 intermediate filament protein as well as the 50 kDa cytokeratin protein designated cytokeratin 14. Basal cells of human epidermis express acidic cytokeratin 14 and basic cytokeratin 5. Cytokeratin 5 is a 58 kDa protein, which is closely related to cytokeratin 6. Cytokeratin 5 and 6 are expressed in basal cell epithelioma, basal cells of prostate, urothelium, vagina, squamous cell carcinomas of skin, tongue, epiglottis and rectal-anal region. Cytokeratin 5 can be useful in the distinction of mesotheliomas from most adenocarcinomas.

Cytokeratin 14 belongs to the type A (acidic) subfamily of low MW cytokeratins and exists in combination with cytokeratin 15. This antibody can be useful to distinguish stratified epithelial cells from simple epithelial cells.

Format

This product is supplied as an antibody cocktail containing sodium azide as a preservative.

Principles of the Procedures

Antigen detection by immunohistochemistry (IHC) is a two-step process involving first, the binding of a primary antibody to the antigen of interest, and second, the detection of bound antibody by a chromogen. The primary antibody may be used in IHC using manual techniques or using Diagnostic BioSystems Automated Montage 360™ Staining System or HighLighter Staining System.

Dilution of Primary Antibody

Diagnostic BioSystems ready to use antibodies have been optimized for use with the recommended Diagnostic BioSystems Detection System and do not require further dilution. Further dilution may result in loss of sensitivity. The user must validate any such change.

Diagnostic BioSystems concentrated antibodies must be diluted in accordance with the staining procedure when used with the recommended Diagnostic BioSystems detection system. Use of any detection methods other than the recommended systems and protocols require validation by the user. Antibody dilutions should be appropriately adjusted and verified according to the detection system used.

Materials Required But Not Provided

Some of the reagents and materials required for IHC are not provided. Pretreatment reagents, detection systems, control reagents and other ancillary reagents are available from Diagnostic BioSystems. Please refer to the Diagnostic BioSystems website at www.dbiosys.com

Storage and Handling

Store at 2-8°C. This antibody is suitable for use until the expiration date when stored at 2-8°C. Do not use product after the expiration date printed on vial. If reagents are stored under conditions other than those specified here, they must be verified by the user. Diluted reagents should be used promptly. Unused portions of antibody preparation should be discarded after one day.

The presence of precipitate or an unusual odor indicates that the antibody is deteriorating and should not be used.

Positive and negative controls should be run simultaneously with all patient specimens. If unexpected staining is observed which cannot be explained by variations in laboratory procedures and a problem with the antibody is suspected, contact Diagnostic BioSystems Technical Support at (925) 484-3350, extension 2 or techsupport@dbiosys.com.

Specimen Collection and Preparation


Tissues fixed in 10% formalin are suitable for use prior to paraffin embedding. Consult references (Kiernan, 1981; Sheehan & Hrapchak, 1980) for further details on specimen preparation.

The user is advised to validate the use of the products with their tissue specimens prepared and handled in accordance with their laboratory practices.

Precautions

This product is a single-use, non-sterile, in vitro diagnostic device.

This antibody contains less than 0.1% sodium azide. Concentrations less than 0.1% are not reportable hazardous materials according to U.S. 29 CFR

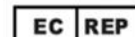

1910.1200, OSHA Hazard Communication and EC Directive 91/155/EC. Sodium azide (NaN₃) used as a preservative is toxic if ingested. Sodium azide may react with lead and copper plumbing to form highly explosive metal azides. Upon disposal, flush with large volumes of water to prevent azide build-up in plumbing. (Center for Disease Control, 1976, National Institute of Occupational Safety and Health, 1976). Specimens, before and after fixation and all materials exposed to them, should be handled as if capable of transmitting infection and disposed of with proper precautions. Never pipette reagents by mouth and avoid contacting the skin and mucous membranes with reagents and specimens. If reagents or specimens come in contact with sensitive areas, wash with copious amounts of water. Microbial contamination of reagents may result in an increase in nonspecific staining. Incubation times or temperatures other than those specified may give erroneous results. The user must validate any such change. The MSDS is available upon request.

Treatment of Tissues Prior to Staining

Place the slides in the recommended antigen retrieval solution using Diagnostic BioSystems Montage Opus™ Antigen Retrieval System if staining is performed manually or with Montage PolyVue Plus™ Auto Detection System. Allow slides to cool down for 20 minutes prior to staining.

Staining Procedure

Refer to the following table for conditions specifically recommended for this antibody. Refer to the Diagnostic BioSystems PolyVue™ Plus—Two Step Detection System or Montage PolyVue Plus™ Auto Detection System or HighLighter core kit for HighLighter Staining System for guidance on specific staining protocols or other requirements.

Parameter	Diagnostic BioSystems Recommendations
Positive Control	Skin, Tonsil
Concentrated Dilution	1:25-1:150
Pretreatment	Citrate Buffer pH 6.0 (Manual/ Montage)
Pretreatment for HL	Tris-EDTA Buffer pH 9.0 (EA1, K102-HL)
Incubation Time & Temperature	30 min @ RT
Detection System	PolyVue™ Plus - Two Step Detection System or Montage PolyVue Plus™ Auto Detection System for Montage 360 System or HighLighter core kit for HighLighter Staining System
Tissue Type	FFPE

Quality Control

Refer to CLSI Quality Standards for Design and Implementation of Immunohistochemistry Assays; Approved Guideline-Second edition (I/LA28-A2) CLSI Wayne, PA USA (www.clsi.org). 2011.

Troubleshooting

Contact Diagnostic BioSystems Technical Support at (925) 484-3350, extension 2, techsupport@dbiosys.com or your local distributor to report unusual staining.

Cellular Localization

Cytoplasmic

Limitations of the Procedure

Immunohistochemistry is a complex technique involving both histological and immunological detection methods. Tissue processing and handling prior to immunostaining can also cause inconsistent results. Variations in fixation and

embedding or the inherent nature of the tissue may cause variations in results (Nadji and Morales, 1983). Endogenous peroxidase activity or pseudoperoxidase activity in erythrocytes and endogenous biotin may cause non-specific staining depending on detection system used. Tissues containing Hepatitis B surface Antigen (HBsAg) may give a false positive with horseradish peroxidase systems (Omata et al, 1980). Improper counterstaining and mounting may compromise the interpretation of results.

Performance Characteristics

The optimum antibody dilution and protocols for a specific application can vary. These include, but are not limited to: fixation, heat-retrieval method, incubation times, and tissue section thickness and detection kit used. Due to the superior sensitivity of these unique reagents, the recommended incubation times and titers listed are not applicable to other detection systems, as results may vary. The data sheet recommendations and protocols are based on exclusive use of Diagnostic BioSystems products. Ultimately, it is the responsibility of the investigator to determine optimal conditions. These products are tools that can be used for interpretation of morphological findings in conjunction with other diagnostic tests and pertinent clinical data by a qualified pathologist.

References

- i) Irvine et al. J Invest Dermatol 108: 809, 1997.
- ii) Whittock et al. Biochem Biophys Res Commun 274: 149, 2000.
- iii) Lane et al. Can Biol 1: 165, 1990.
- iv) Moll et al. Cell 31: 11, 1982.
- v) Purkis et al. J Cell Sci 97: 39, 1990.
- vi) Westzels et al. Am J Pathol 134: 571, 1989

